# REUSE OR EXTENDED USE OF EYE PROTECTION

COVID-19

### **REUSE:**

removal and storage of eye protection between patients



### **EXTENDED USE:**

use of the same eye protection for repeated close contact with multiple patients without removal of eye protection between patients (acceptable when patients have the same respiratory pathogen)


# Considerations for extended use or reuse of <u>eye protection</u> to reduce contact transmission:

Acceptable for <u>disposable or reusable devices</u>

 If eye protection becomes soiled or difficult to see through, it should be removed and reprocessed; discard if damaged

 Leave patient care area if eye protection must be removed

 Do not touch eye protection--if eye protection is touched, wearer should immediately perform proper hand hygiene measures

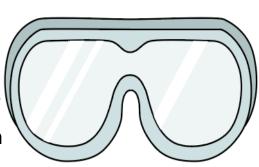

## Considerations for reprocessing eye protection:

- Adhere to manufacturer instructions if available
- If instructions are unavailable, use the following guidelines:
  - (1) While wearing gloves, wipe the inside, followed by the outside using a clean cloth with neutral detergent or a cleaner wipe
  - (2) Wipe the outside with a wipe or clean cloth saturated in EPAregistered hospital disinfectant solution
  - (3) Wipe outside with clean water or alcohol, removing residue
  - (4) Fully dry (air dry or with clean, absorbent towels)
  - (5) Remove gloves from hands and perform hand hygiene

### **Storage & Donning of Used Eye Protection:**

- Eye protection should be stored in a clean sealable paper bag or other breathable container
- Reused or extended use eye protection should have ONLY ONE WEARER
- Clean hands with soap and water OR alcohol-based sanitizer before and after touching/adjusting the eye protection for wear
- Do not touch the outside of the eye protection--if contact occurs, wash hands as described above
- <u>Wear a pair of non-sterile gloves when re-donning</u> the eye protection, discard gloves when finished



#### For more information

visit our website at www.gchd.us/coronavirus and follow us on Facebook!

Last reviewed: 06/08/2020